


Simplicity	inch	mm
A Minimum swivel radius	26¾	680
B Minimum wall to seat back	3	75
C Seat depth	141/2	370
D Minimum folded width to footplate	111/4	285
E Minimum open width to edge of footplate	22¾	580
F Armrest width – external	241/2	625
G Armrest width - internal	19	480
H Seat back height	15¾	400
Footplate to seat height standard seat option extended seat post	15¼-19¼ 18-22	390-490 460-560
J Minimum footplate height	2	50
Bottom parking space required	26	660 [†]
Minimum track intrusion into staircase	61/4	160
Minimum track intrusion into top landing (typical at 45°)	6	150*
Maximum weight limit – hinge/slidetrack	309lb / 140.45 kgs 265lb / 120.45 kgs	
Speed	4%"m/sec (approx)	
Power supply	24v DC battery	
Charger	240v	
Staircase angle	27-50°	

†Increases by up to 140mm with zero intrusion kit.

*A zero track intrusion kit is available.



Sterling Stairlifts from Handicare are backed by a 2 year parts warranty plus a LIFETIME warranty on the gearbox and motor.

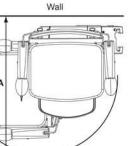

Sterling Stairlifts - Handicare Accessibility 10 Industrial Park Drive

Mount Pocono, PA 18344 Toll Free: 866-276-5438 Local: 570-895-2480 Fax: 570-895-2481

info@sterling-stairlifts.com www.sterling-stairlifts.com



Please be aware that actual colors may vary slightly from those shown. All information is provided in good faith but does not constitute a guarantee. Handicare reserves the right to alter designs and specifications without prior notice.







SimplicityStraight stairlift



1 E. Hollis St, Nashua, NH 03060 800-660-7221 www.nhlift.com

The benefits of a Sterling Simplicity Stairlift from Charron Medical Services

If you are thinking about installing a stairlift the Simplicity offers you a safe and reliable solution.

Handicare's Sterling Stairlifts have helped thousands of people overcome the problems posed by going up and down their stairs at home.

The Simplicity straight stairlift offers you all the basic functions that you would expect from a stairlift and, if you have a doorway close to the base of the stairs, you can choose from a manually operated hinged track or the unique SlideTrack

that is only available from Handicare.

Supplied in a neutral shade, the wipe clean upholstery will

compliment most color schemes.

The track is fitted to the stairs, not the wall, so you will not need to redecorate following the installation of your Simplicity stairlift.

Like all Handicare Sterling stairlifts the Simplicity has a battery backup which means it will continue to work in the event of a power outage so you do not need to worry about being stranded up or downstairs.

The compact unit folds up leaving the staircase free for other users and also comes with two remote controls which means, once folded, you can move it along the track freeing up space in your hallway. The Simplicity will charge its batteries anywhere on the track that you choose to leave it.









Safety As Standard

Your Simplicity stairlift will come with a range of safety sensors that are designed to stop your stairlift if it meets any obstruction as you travel up the stairs.

Your Simplicity stairlift has a clip operated seat belt for added security and peace of mind.

Smooth & Quiet

Because of the high quality motors we use in our Simplicity stairlifts we are able to offer perhaps the quietest and smoothest stairlift ride available. All backed by a lifetime gearbox and motor guarantee.

Easy to Use

The operating toggle is easy to use, simply hold it towards the direction you wish to travel, up or down, and your Simplicity will move smoothly and gently on its way.

The control is conveniently located at the end of the arm leaving both arms clear to provide unhindered support when getting on or off the seat.

